REAL TIME PCR DETECTION KIT

Campylobacter

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VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit 6 x 8-well strips, low profile	VS-CAM106L
VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit 6 x 8-well strips, high profile	VS-CAM106H
VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit 12 x 8-well strips, low profile	VS-CAM112L
VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit 12 x 8-well strips, high profile	VS-CAM112H
VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit 96-well plate, low profile	VS-CAM113L
VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit 96-well plate, high profile	VS-CAM113H

VIASURE

ENGLISH

1. Intended of use

VIASURE *Campylobacter* Real Time PCR Detection Kit is designed for specific identification of *Campylobacter* in human stool samples from patients with signs and symptoms of gastrointestinal infection. This test is intended for use as an aid in the diagnosis of *Campylobacter* infection in humans in combination with clinical and epidemiological risk factors. DNA is extracted from stool specimens, multiplied using Real Time amplification and detected using specific primers and a fluorescent reporter dye probe for *Campylobacter*.

2. Summary and Explanation

The *Campylobacter* genus belongs to the family *Campylobacteraceae*. Up to date, it consists of 26 species, 2 provisional species, and 9 subspecies. The most common species of *Campylobacter* associated with human illness are *C. jejuni* and *C. coli*, but other species can also cause human infections.

Campylobacteriosis is the most prevalent foodborne bacterial infection increased in both developed and developing countries. Poultry is a major reservoir and source of transmission of campylobacteriosis to humans. Other risk factors include consumption of animal products and water, contact with animals, and even person-to-person transmission (fecal-oral or via fomites).

Infection with *Campylobacter* causes gastroenteritis characterised by fever, vomiting, headaches, and abdominal pain with watery or bloody diarrhea, for a median duration of 6 days. However, in a minority of individuals *Campylobacter* infection may lead to autoimmune conditions known as Guillain-Barré syndrome (GBS) and Miller Fisher syndrome (MFS). In addition, they could be involved in extragastrointestinal manifestations, including bacteremia, lung infections, brain abscesses, meningitis, and reactive arthritis.

The Real Time PCR assays have used extensively the 16S rRNA gene for rapid detection and identification of many bacterial taxa, including *Campylobacter* species.

3. Principle of the procedure

VIASURE *Campylobacter* Real Time PCR Detection Kit is designed for the diagnosis of gastroenteritis caused by *Campylobacter* in human stool samples. After DNA isolation, the identification of *Campylobacter* is performed by the amplification of a conserved region with the 16S rRNA gene, using specific primers and a fluorescent-labeled probe.

VIASURE *Campylobacter* Real Time PCR Detection Kit is based on 5' exonuclease activity of DNA polymerase. During DNA amplification, this enzyme cleaves the probe bound to the complementary DNA sequence, separating the quencher dye from the reporter. This reaction generates an increase in the fluorescent signal which is proportional to the quantity of the hydrolyzed target sequence. This fluorescence could be measured on real time PCR platforms.

VIASURE *Campylobacter* Real Time PCR Detection Kit contains in each well all the components necessary for real time PCR assay (specific primers/probes, dNTPS, buffer, polymerase) in a stabilized format, as well as an internal control to discard the inhibition of the polymerase activity. *Campylobacter* DNA targets are amplified and detected in FAM channel and the internal control (IC) in HEX, VIC or JOE channel (depending on the equipment used).

4. Reagents provided

VIASURE *Campylobacter* Real Time PCR Detection Kit includes the following materials and reagents detailed in Table 1 and Table 2:

Reference	Reagent/Material	Description	Color	Amount
VS-CAM1SL/ VS-CAM1SH	<i>Campylobacter</i> 8-well strips	A mix of enzymes, primers probes, buffer, dNTPs, stabilizers and Internal control in stabilized format	White	6/12 X 8-well strip
VS-RB02	Rehydration Buffer	Solution to reconstitute the stabilized product	Blue	1 vial x 1.8 mL
VS-CAM1C	<i>Campylobacter</i> Positive Control	Non-infectious synthetic lyophilized DNA	Red	1 vial
VS-NC1	Negative control	Non template control	Violet	1 vial x 1 mL
VS-H2O	Water RNAse/DNAse free	Water RNAse/DNAse free	White	1 vial x 1 mL
VS-OCS	Tear-off 8-cap strips	Optical caps for sealing wells during thermal cycling	Transparent	6/12 X 8-cap strip

Table 1. Reagents and materials provided in VIASURE *Campylobacter* Real Time PCR Detection Kit with Ref. VS-CAM106L, VS-CAM106H, VS-CAM112L and VS-CAM112H.

Reference	Reagent/Material	Description	Color	Amount
VS-CAM1PL/ VS-CAM1PH	<i>Campylobacter</i> 96-well plate	A mix of enzymes, primers probes, buffer, dNTPs, stabilizers and Internal control in stabilized format	White	1 plate
VS-RB02	Rehydration Buffer	Solution to reconstitute the stabilized product	Blue	1 vial x 1.8 mL
VS-CAM1C	<i>Campylobacter</i> Positive Control	Non-infectious synthetic lyophilized DNA	Red	1 vial
VS-NC1	Negative control	Non template control	Violet	1 vial x 1 mL
VS-H2O	Water RNAse/DNAse free	Water RNAse/DNAse free	White	1 vial x 1 mL
VS-OCS	Tear-off 8-cap strips	Optical caps for sealing plate during thermal cycling	Transparent	12 X 8-cap strip

Table 2. Reagents and materials provided in VIASURE *Campylobacter* Real Time PCR Detection Kit with Ref. VS-CAM113L and VS-CAM113H.

5. Reagents and equipment to be supplied by the user

The following list includes materials that are required for using, but not included in the VIASURE *Campylobacter* Real Time PCR Detection Kit.

- Real Time PCR instrument (thermocycler) (to check compatibility see Annex I).
- DNA extraction kit.
- Centrifuge for 1.5 mL tubes.
- Vortex.
- Micropipettes (0.5-20 µL, 20-200 µL).
- Filter tips.
- Powder-free disposal gloves.

6. Transport and storage conditions

- The kits can be shipped at 2-50°C for a maximum of 2 weeks.
- For long term storage the kit should be stored at 2-8°C.
- Once the positive control has been re-suspended, store it at -20°C. We recommend to separate in aliquots to minimize freeze and thaw cycles.
- Keep components away from sunlight.

7. Precautions for users

- For professional *in vitro* diagnostic use.
- Do not use after expiration date.
- Design a unidirectional process flow. It should begin in the Extraction Area and then move to the Amplification and Detection Area. Do not return samples, equipment and reagents to the area in which the previous step was performed.
- Follow Good Laboratory Practices. Wear protective clothing, use disposal gloves, goggles and mask. Do not eat, drink or smoke in the working area. Once you finish the test wash your hands.
- Specimens must be treated as potentially infectious as well as all reagents and materials that have been exposed to the samples and must be handled according to the national safety regulations. Take necessary precautions during the collection, storage, treatment and disposal of samples.
- Regular decontamination of commonly used equipment is recommended, especially micropipettes and work surfaces.

8. Test procedure

8.1. SAMPLE PREPARATION

Stool samples should be collected in clean containers and processed as soon as possible to guarantee the quality of the test. We recommend to use fresh samples.

For longer storage, the samples must be frozen at -20°C. In this case, the sample will be totally thawed and brought to room temperature before testing. Homogenise stool sample as thoroughly as possible prior to preparation. Freezing and thawing cycles are not recommended.

We recommend to dilute stool samples before extraction. Collect a pea-size stool (approx. 8mm) and put in a 1.5 mL microcentrifuge tube containing 100 µL of physiological serum or PBS. Vortex intensely and centrifuge 10,000 rpm for 1min. Use 200 µL of supernatant to perform DNA extraction.

8.1.1. DNA EXTRACTION

For DNA extraction from human stool samples you can use your manually or automatic routine optimized system. Also, you can use any commercially available DNA extraction kit and follow the manufacturer's instructions for use. We have validated the following extraction kits:

- Viasure RNA-DNA Extraction kit (VIASURE), recommended.
- Invisorb® Spin Universal Kit (Stratec).
- UltraClean® Tissue & Cells DNA/RNA Isolation Kit (Mobio).
- NucleoSpin® RNA Virus (Machery Nagel).

8.2. LYOPHILIZED POSITIVE CONTROL

Campylobacter Positive Control contains high copies template, the recommendation is to open and manipulate it in a separate laboratory area away from the other components. Reconstitute the lyophilized *Campylobacter* Positive Control (red vial) adding 100 µL of Water RNase/DNase free (white vial) supplied and vortex thoroughly.

Once the positive control has been re-suspended, store it at -20°C. Recommendation is to separate it in aliquots to minimize freeze and thaw cycles.

8.3. PCR PROTOCOL

Determine and separate the number of required reactions including samples and controls. One positive and negative control must be included in each run. Peel off protective aluminum seal from plates or strips.

1) Reconstitute the number of wells you need.

Add 15 µL of Rehydration Buffer (blue vial) into each well.

2) Adding samples and controls.

Add 5 µL of DNA sample, reconstituted *Campylobacter* Positive Control (red vial) or Negative Control (violet vial) in different wells and close the wells with the caps provided. Centrifuge briefly.

Load the plate or the strips in the thermocycler.

3) Set up your thermocycler.

Program your thermocycler following the conditions below and start the run:

Cycles	Step	Time	Temperature
1	Polymerase activation	2 min	95°C
45	Denaturalization	10 seg	95°C
	Annealing/Data collection*	50 seg	60°C

Table 3. PCR protocol

Fluorogenic data should be collected during the annealing step (*) through the FAM (*Campylobacter*) and HEX, JOE or VIC channels (Internal Control (IC)). In Applied Biosystems 7500 Fast Real-Time PCR System, Applied Biosystems StepOne™ Real-Time PCR System and Stratagene Mx3005P™ Real Time PCR System check that passive reference option ROX is none.

9. Result interpretation

The use of positive and negative controls in each run, validates the reaction by checking the absence of signal in negative control well and the presence of signal for *Campylobacter* positive control well. Check Internal Control signal to verify the correct functioning of the amplification mix. The analysis of the samples is done by the software itself of the used real time PCR equipment according to manufacturer’s instructions for use.

Using the following table, read and analyze the results:

<i>Campylobacter</i>	Internal control	Negative control	Positive control	Interpretation
+	+/-	-	+	<i>Campylobacter</i> Positive
-	+	-	+	<i>Campylobacter</i> Negative
+	+	+	+	Experiment fail
-	-	-	-	Experiment fail

Table 4. Sample interpretation


+: Amplification curve
 -: No amplification curve

A sample is considered positive if the Ct value obtained is less than 40 and the internal control shows an amplification signal.

A sample is considered positive if the sample shows an amplification signal less than 40 Ct value but the internal control is negative. Sometimes, the detection of internal control is not necessary because a high copy number of target can cause preferential amplification of target-specific nucleic acids.

A sample is considered negative, if the sample shows no amplification signal in the detection system but the internal control is positive. An inhibition of the PCR reaction can be excluded by the amplification of internal control.

Figure 1. Correct run of negative and positive control run on the Bio-Rad CFX96 Touch™ Real-Time PCR Detection System.

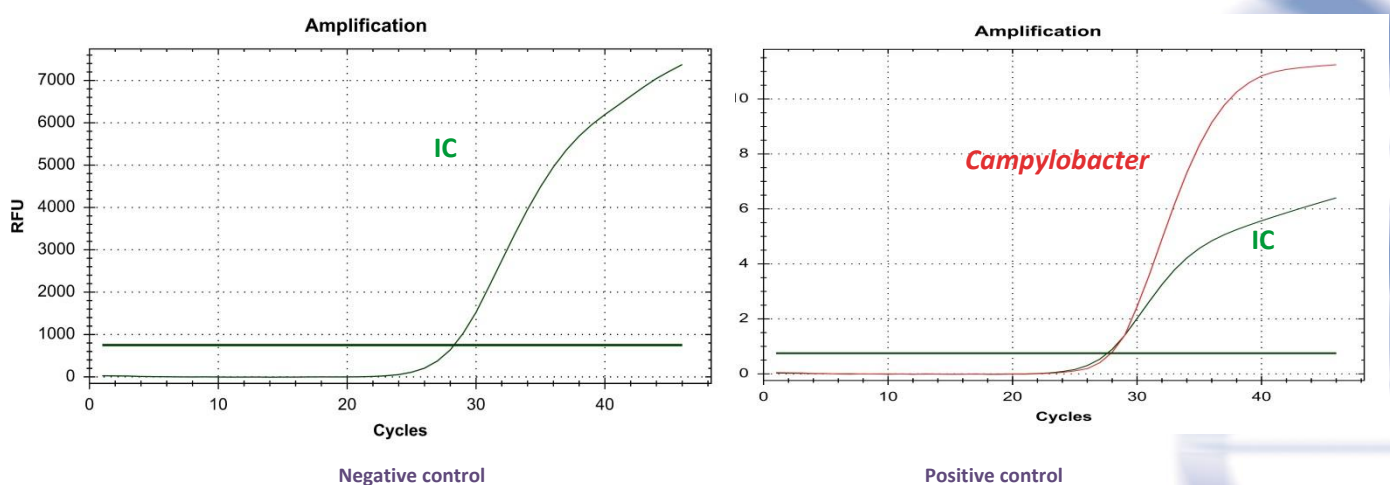The result is considered invalid if there is signal of amplification in negative control or absence of signal in the positive well. We recommend to repeat the assay again.

In case of absence of internal control signal in sample wells we recommend to repeat the assay diluting the sample 1:10 or to repeat the extraction to check for possible problems of inhibition.

10. Limitations of the test

- The result of the test should be evaluated by a health care professional and evaluated in the context of medical history, clinical symptoms and other diagnostic tests.
- Although this assay can be used with other types of samples it has been validated with human faecal samples.
- The quality of the test depends on the quality of the sample; proper DNA from faecal specimens must be extracted. Unsuitable collection, storage and/or transport of specimens may give false negative results.
- Extremely low levels of target below the limit of detection may be detected, but results may not be reproducible.
- There is a possibility of false positive results due to cross-contamination by *Campylobacter*, either samples containing high concentrations of target DNA or contamination due to PCR products from previous reactions.

11. Quality control

VIASURE *Campylobacter* Real Time PCR Detection Kit contains a positive and a negative control that must be included in each run to correctly interpret the results. Also, the internal control (IC) in each well confirms the correct performance of the technique.

12. Performance characteristics

12.1. CLINICAL SENSITIVITY AND SPECIFICITY

Evaluation with faecal samples from symptomatic patients was tested using VIASURE *Campylobacter* Real Time PCR Detection Kit. These results were compared with those obtained by a commercial Real Time PCR Kit (RIDA®GENE Bacterial Stool Panel (r-Biopharm)).

The results were as follows:

VIASURE <i>Campylobacter</i> Real Time PCR Detection Kit	RIDA®GENE Bacterial Stool Panel (r-Biopharm)			Total
		+	-	
+	42	3 ^{*#}	45	
-	0	59	59	
Total	42	62	104	

Table 5. Comparative results

* The low amount of template DNA detected in these faecal samples is below the detection limit of the method used and so it may produce random positive results.


Besides, these two samples have been evaluated by a conventional PCR and subsequently sequencing, confirming our results.

The results show a high sensitivity and specificity to detect *Campylobacter* using VIASURE *Campylobacter* Real Time PCR Detection Kit.

12.2. ANALYTICAL SENSITIVITY

VIASURE *Campylobacter* Real Time PCR Detection Kit has a detection limit of ≥ 10 DNA copies per reaction (Figure 2).

Figure 2. Dilution series of *Campylobacter* (10^7 - 10^1 copies/rxn) template run on the Bio-Rad CFX96 Touch™ Real-Time PCR Detection System.

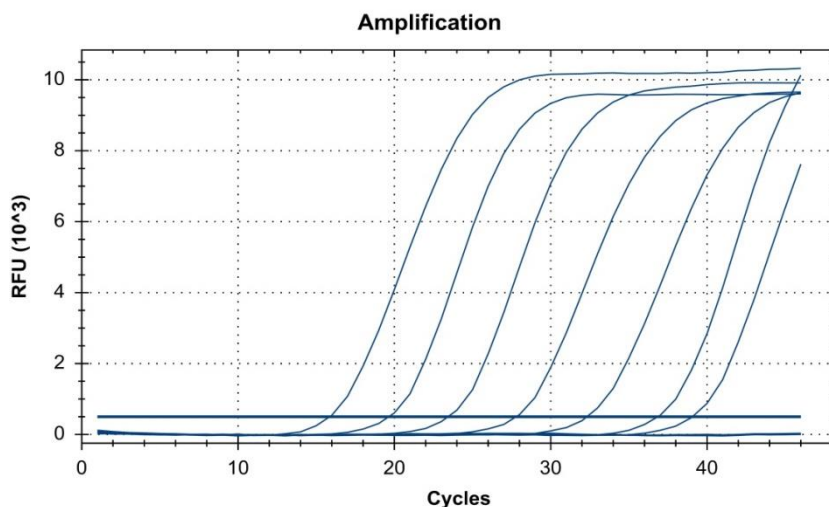12.3. ANALYTICAL SPECIFICITY

The specificity of the *Campylobacter* assay was confirmed by testing a panel consisting of 40 microorganisms representing the most common enteric pathogens or flora present in the intestine.

VIASURE *Campylobacter* Real Time PCR Detection Kit not identify the following strains confirming the exclusivity of the trial.

Cross-reactivity testing					
<i>Shigella flexneri</i>	-	<i>Helicobacter heilmannii</i>	-	<i>Arcobacter butzleri</i>	-
<i>Shigella dysenteriae</i>	-	<i>Proteus vulgaris</i>	-	<i>Pseudomonas aeruginosa</i>	-
<i>Salmonella typhi</i>	-	<i>Aeromonas hydrophila subsp. hydrophila</i>	-	<i>Enterococcus faecalis</i>	-
<i>Salmonella paratyphi A</i>	-	<i>Citrobacter freundii</i>	-	<i>Giardia intestinalis</i>	-
<i>Salmonella paratyphi B</i>	-	<i>Staphylococcus aureus subsp. aureus</i>	-	<i>Yersinia enterocolitica O:3</i>	-
<i>Salmonella typhimurium</i>	-	<i>Serratia liquefaciens</i>	-	<i>Yersinia enterocolitica O:9</i>	-
<i>Salmonella bongori</i>	-	<i>Vibrio parahaemolyticus</i>	-	<i>Bacteroides fragilis</i>	-
<i>Salmonella enteritidis</i>	-	<i>Clostridium difficile</i>	-	<i>Cryptosporidium parvum</i>	-
<i>Salmonella enterica subsp. entérica</i>	-	<i>Clostridium perfringens</i>	-	<i>Entamoeba histolytica</i>	-
<i>Salmonella pullorum</i>	-	<i>Enterotoxigenic Escherichia coli</i>	-	Adenovirus serotype 40	-
<i>Salmonella gallinarum</i>	-	<i>Enteropathogenic Escherichia coli</i>	-	Adenovirus serotype 41	-
<i>Helicobacter pylori</i>	-	<i>Klebsiella oxytoca</i>	-	Rotavirus A	-
<i>Helicobacter hepaticus</i>	-	<i>Listeria monocytogenes</i>	-	Norovirus Genotypes I and II	-
<i>Helicobacter cinaedi</i>	-	<i>Candida albicans</i>	-	Astrovirus Genotype I-VIII	-

Table 6. Reference pathogenic and non-pathogenic microorganisms used in this study.

12.4. ANALYTICAL REACTIVITY

The reactivity of VIASURE *Campylobacter* Real Time PCR Detection Kit was evaluated against *Campylobacter jejuni subsp. jejuni*, *Campylobacter coli*, *Campylobacter lari*, *Campylobacter upsaliensis* and *Campylobacter fetus* showing positive results.

ANNEX 1:

COMPATIBILITY OF THE MOST COMMON REAL TIME PCR EQUIPMENT

Low profile strips can be used in all PCR thermocyclers equipped with low profile block, like systems listed in table A.1. High profile strips can be used in all PCR thermocyclers equipped with high or regular profile block, like systems listed in table A.2. If you do not find your thermocycler in the list below, please contact with your supplier.

Table A.1 LOW PROFILE BLOCK THERMOCYCLERS	
Manufacturer	Model
Applied Biosystems	7500 Fast Real-Time PCR System
Applied Biosystems	7500 Fast Dx Real-Time PCR System
Applied Biosystems	QuantStudio™ 12K Flex 96-well Fast
Applied Biosystems	QuantStudio™ 6 Flex 96-well Fast
Applied Biosystems	QuantStudio™ 7 Flex 96-well Fast
Applied Biosystems	StepOne Plus™ Real-Time PCR System
Applied Biosystems	StepOne™ Real-Time PCR System
Applied Biosystems	ViiA™ 7 Fast Real-Time PCR System
Bio-Rad	CFX96 Touch™ Real-Time PCR Detection System
Bio-Rad	Mini Opticon™ Real-Time PCR Detection System
Roche	LightCycler®480 Real-Time PCR System
Roche	LightCycler®96 Real-Time PCR System
Agilent Technologies	AriaMx Realtime PCR System

Table A.2 HIGH PROFILE BLOCK THERMOCYCLERS	
Manufacturer	Model
Applied Biosystems	7300 Real-Time PCR System
Applied Biosystems	7500 Real-Time PCR System
Applied Biosystems	7900 HT Real-Time PCR System
Applied Biosystems	ABI PRISM 7000
Applied Biosystems	ABI PRISM 7700
Applied Biosystems	QuantStudio™ 12K Flex 96-well
Applied Biosystems	QuantStudio™ 6 Flex 96-well
Applied Biosystems	QuantStudio™ 7 Flex 96-well
Applied Biosystems	ViiA™ 7 Real-Time PCR System
Bio-Rad	CFX96 Touch™ Deep Well Real-Time PCR Detection System
Bio-Rad	iCycler iQ™ Real-Time PCR Detection System
Bio-Rad	iCycler iQ™5 Real-Time PCR Detection System
Bio-Rad	MyiQ™ Real-Time PCR Detection System
Bio-Rad	MyiQ™2 Real-Time PCR Detection System
Eppendorf	Mastercycler™ep <i>realplex</i>
Stratagene / Agilent Technologies	Mx3000P™ Real Time PCR System
Stratagene / Agilent Technologies	Mx3005P™ Real Time PCR System
Analytik Jena Biometra	TOptical
Analytik Jena Biometra	qTOWER 2.0

Table A1/A2. Compatible low and high profile Real Time PCR systems.

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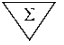
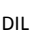
13. Bibliography/Bibliografía

1. N.O. Kaakoush et al. Global Epidemiology of *Campylobacter* Infection. *Clinical Microbiology Reviews* 2015; 28(3): 687-720.
2. S.M. Man. The clinical importance of emerging *Campylobacter* species. *Nature Reviews Gastroenterology & Hepatology* 2011; 8(12): 669-685.
3. R.F. de Boer et al. Detection of *Campylobacter* species and *Arcobacter butzleri* in stool samples by use of real-time multiplex PCR. *Journal of Clinical Microbiology* 2013; 51(1): 253-259.

14. Symbols for IVD components and reagents/Símbolos para reactivos y productos para diagnóstico *in vitro*

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
VIASURE *Campylobacter* Real Time PCR Detection Kit has been validated on the following equipments: Applied Biosystems 7500 Fast Real-Time PCR System, Applied Biosystems StepOne™ Real-Time PCR System, Bio-Rad CFX96 Touch™ Real-Time PCR Detection System and AriaMx Realtime PCR System. When using the Applied Biosystems 7500 Fast with strips it is recommend to place a plate holder to reduce the risk of crushed tube (Ref. PN 4388506).

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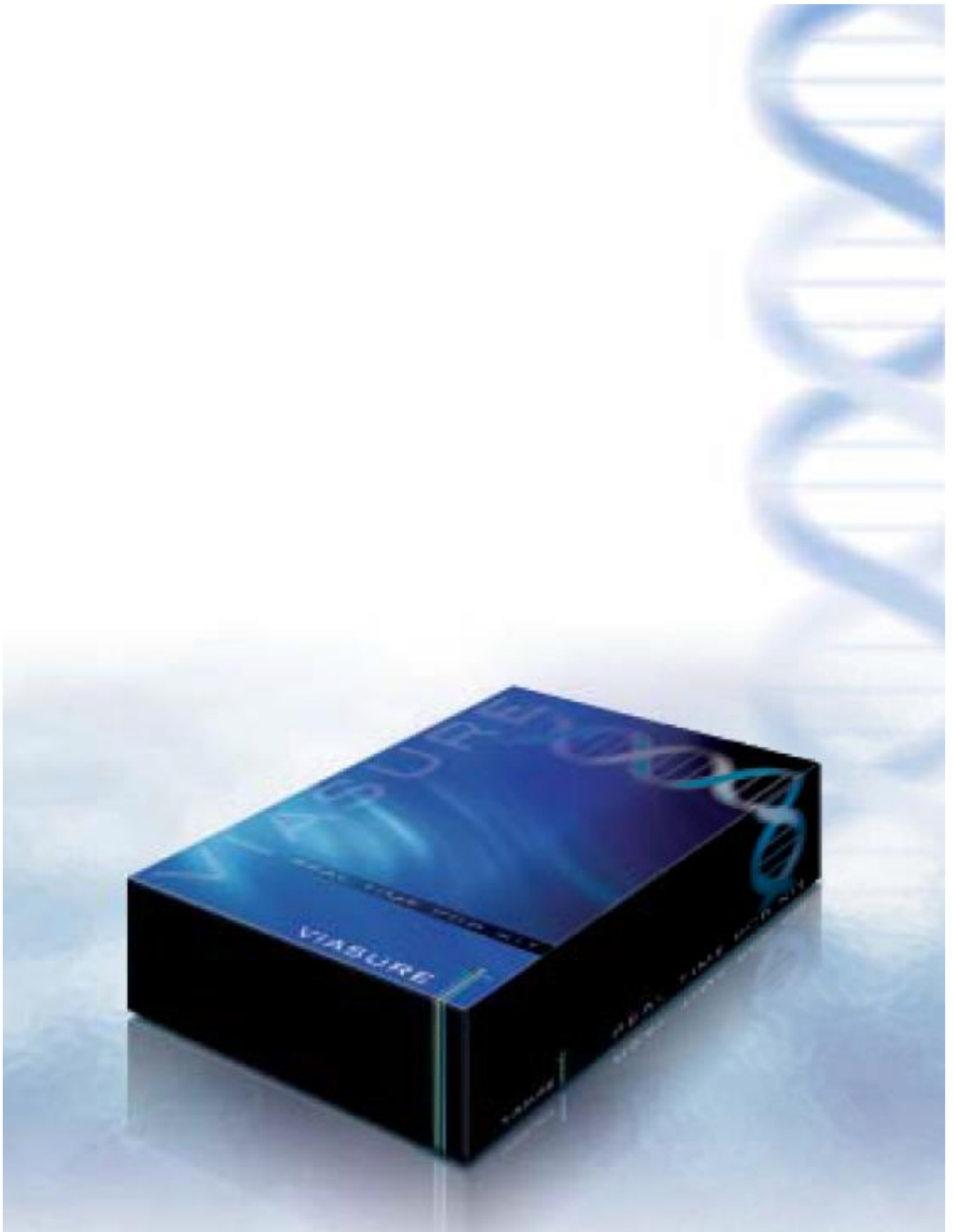- CFX™ and IQ5™ are registered trademarks of Bio-Rad Laboratories.
- ABI®, QuantStudio™, StepOnePlus™ and ViiA™ are registered trademarks of Thermo Fisher Scientific Inc.
- LightCycler® is a registered trademark of Roche.
- Mx3000P™, Mx3005™ and AriaMx Realtime PCR System are registered trademarks of Agilent Technologies.
- Mastercycler™ is a registered trademark of Eppendorf.

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